




User Manual

Rouge & Rouge Pro Series Therapy Lights



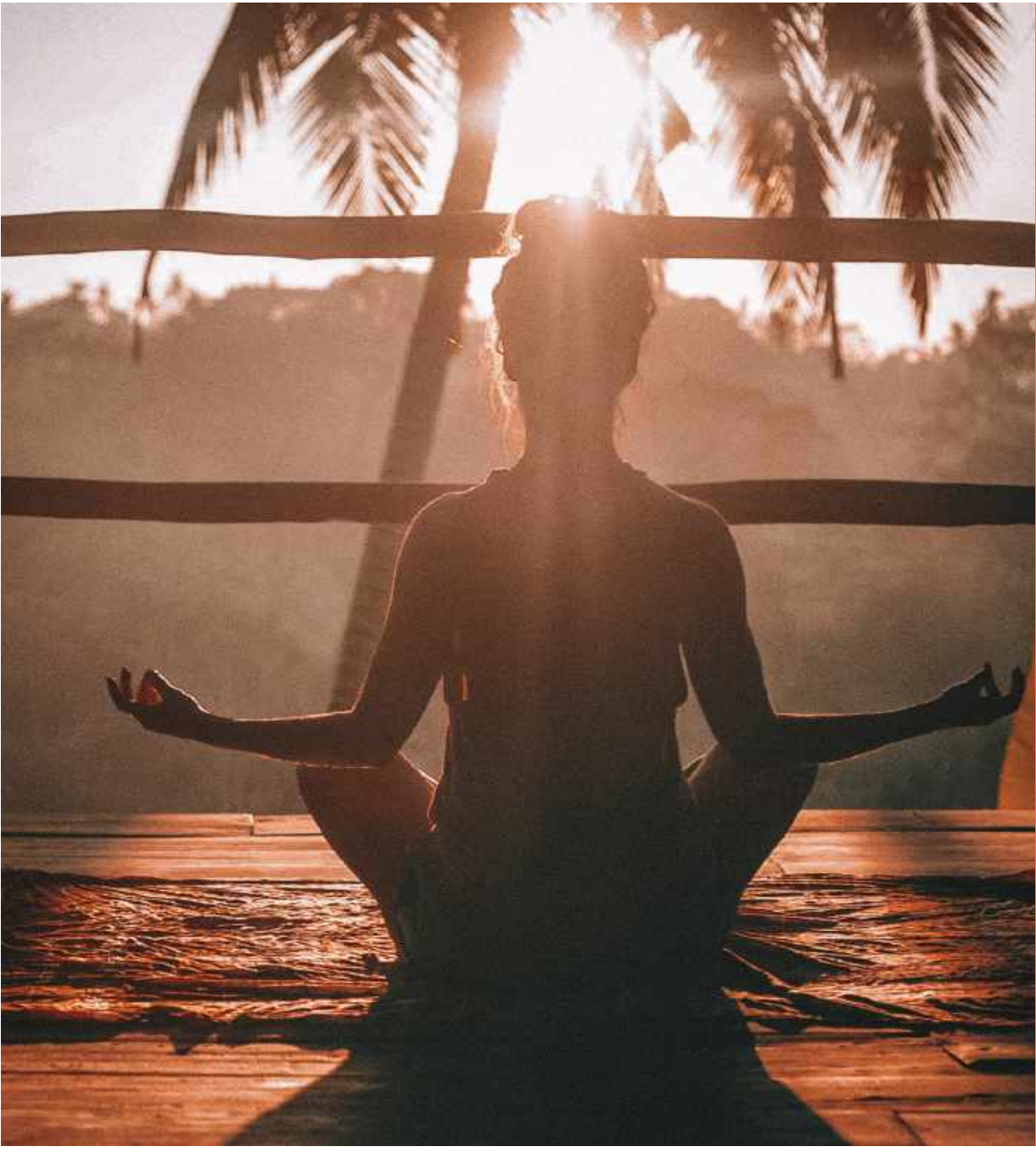

Healthy Lifestyle Choice



TABLE OF CONTENTS


Understanding Red Light Therapy	/01	Why choose us	/02
Rouge Pro Series Compare Specs	/04	Rouge Pro Series Specifications	/05
Rouge Series Compare Specs	/07	Rouge Series Specifications	/08
Rouge Pro Series Irradiance	09/10	Rouge Pro Series Optional Mounting Configurations and Instructions	/11
Rouge Series Optional Mounting Configurations and Instructions	/12	Package Contents	/13
Install Your Therapy Light	14/15	Setting up Your Rouge Pro Series Light	/16
Connecting Multiple Rouge Pro Series Light Units	17/18	Setting up Your Rouge Series Light	/19
Connecting Multiple Rouge Series Light Units	/20	Warranty	/21

What is Red Therapy Light and How Does it Work?

Red light therapy is also known as low-level laser therapy (LLLT). It uses low-energy light emitting diodes (LEDs) or lasers as an alternative therapy for pain relief, or to promote good cell functions. The principle of red light therapy can be thought of as like a human version of the way plants photosynthesize. It is thought that certain

wavelengths of light can enhance the functioning of mitochondria, which are the organelles in a cell which take in nutrients and break them down to create energy molecules for the cell. By enhancing the cellular respiration, it will improve the functioning of the whole cell.


Benefits Of Rouge LED Light Therapy

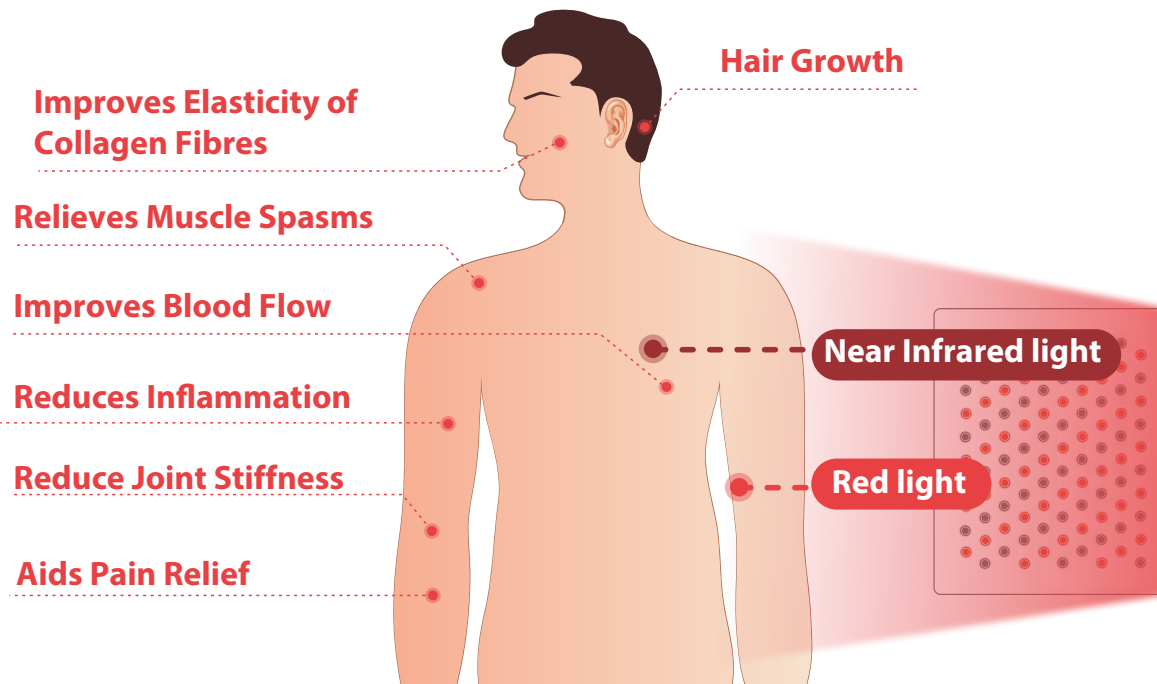



Why choose us?

1. **More Power:** Rouge panels are the most powerful home use devices on the market today. Compare us to the competition and you'll notice we deliver up to 5 times more power per comparable panel.
2. **More value:** our products are priced for home use, while our competitors price their panels for medical clinics and sports professionals.
3. **Satisfaction Guaranteed:** If for any reason you're not totally satisfied with your panel, you can return it to us within 60 days for a full refund.

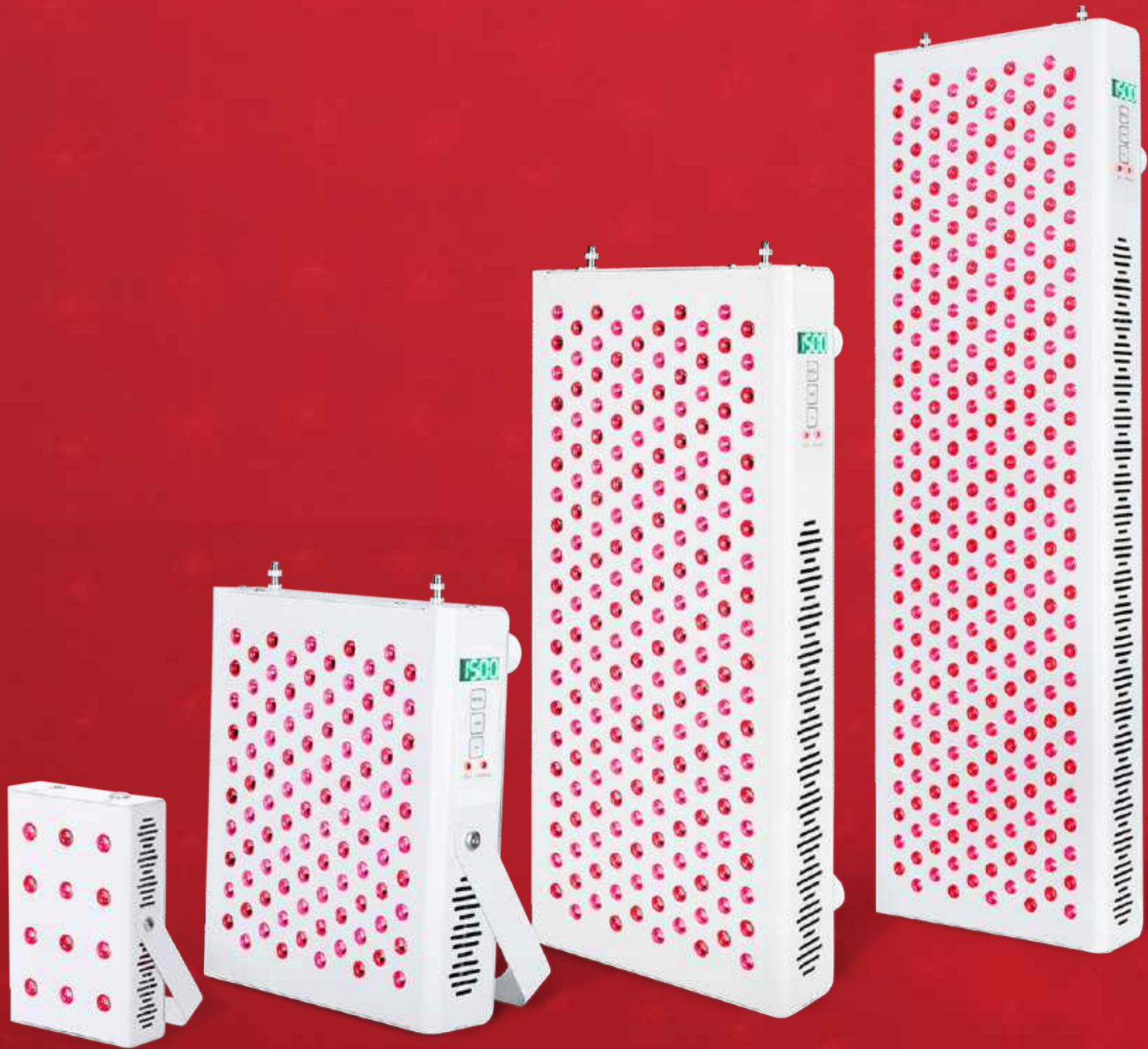

Rouge Pro Series Product Specifications

	TRAVEL EDITION	PRO 500	PRO 1000	PRO 1500
LED Power	60W	500W	1000W	1500W
Wavelengths:	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm
Dimensions	5.9" x 3.9" x 1.4"	10.4" x 8.7" x 2.3"	18.9" x 8.7" x 2.3"	27.2" x 8.7" x 2.3"
Max Area	—	36" x 34"	46" x 34"	52" x 34"
Weight	0.85KG/1.9LBS	2.8KG/6.2LBS	4.5KG/9.9LBS	7KG/15.4LBS
Irradiance (mW/cm2)*	> 100 mw/cm2 at the surface of the LEDs	135	145	161
EMF Emission	– Negligible	– Negligible	– Negligible	– Negligible
Warranty	3 Years	3 Years	3 Years	3 Years
Daisy-chained	<input type="checkbox"/>	Maximum 6pcs	Maximum 5pcs	Maximum 4pcs
Wireless & Rechargeable	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Door Mount	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Mobile Stand	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Tabletop Stand	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
WiFi Ready	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Built-in Programmable Timer	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Smart Built-in Cooling Fans	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>

Rouge Pro Series Product Specifications

	TRAVEL EDITION	PRO 500	PRO 1000	PRO 1500
LED Power	60W	500W	1000W	1500W
Wavelengths	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm
Amps @110v	—	1.1A	2.1A	3.3A
Amps @220v	—	0.4A	1.1A	1.3A
Dimensions	5.9" x 3.9" x 1.4"	10.4" x 8.7" x 2.3"	18.9" x 8.7" x 2.3"	27.2" x 8.7" x 2.3"
Power Consumption	14W	120W	235W	352W
LEDs	12 x 5W	100 x 5W	200 x 5W	300 x 5W
Beam Angle	30 Degree	30 Degree	30 Degree	30 Degree
Cooling Fans	No	1pc fan	3pcs fans	4pcs fans
Weight	0.85KG/1.9LBS	2.8KG/6.2LBS	4.5KG/9.9LBS	7KG/15.4LBS
Lifespan	100,000hrs	100,000hrs	100,000hrs	100,000hrs
EMF Emission	– Negligible	– Negligible	– Negligible	– Negligible
Warranty	3 Years	3 Years	3 Years	3 Years
Best For	Targeted Treatment	Targeted Treatment	Targeted Treatment	Full Body Treatment
Max Area	—	36" x 34"	46" x 34"	52" x 34"





Rouge Series Product Specifications

	ROUGE 500	ROUGE 1000	ROUGE 1500
LED Power	500W	1000W	1500W
Wavelengths	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm
Dimensions	10.4" x 8.7" x 2.3"	18.9" x 8.7" x 2.3"	27.2" x 8.7" x 2.3"
Max Area	36" x 34"	46" x 34"	52" x 34"
Weight	2.8KG/6.2LBS	4.5KG/9.9LBS	7KG/15.4LBS
Irradiance (mW/cm2)*	135	145	161
EMF Emission	-Negligible	-Negligible	-Negligible
Warranty	3 Years	3 Years	3 Years
Daisy-chained	Maximum 6pcs	Maximum 5pcs	Maximum 4pcs
Wireless & Rechargeable	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Door Mount	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Mobile Stand	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Tabletop Stand	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
WiFi Ready	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Built-in Programmable Timer	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Smart Built-in Cooling Fans	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>


Rouge Series Product Specifications

ROUGE 500

ROUGE 1000


ROUGE 1500

LED Power	500W	1000W	1500W
Wavelengths	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm	42mil RED 660nm & NIR 850nm
Amps @110v	1.1A	2.1A	3.3A
Amps @220v	0.4A	1.1A	1.3A
Dimensions	10.4" x 8.7" x 2.3"	18.9" x 8.7" x 2.3"	27.2" x 8.7" x 2.3"
Power Consumption	120W	235W	352W
LEDs	100 x 5W	200 x 5W	300 x 5W
Beam Angle	30 Degree	30 Degree	30 Degree
Cooling Fans	1pc fan	3pcs fans	4pcs fans
Weight	2.8KG/6.2LBS	4.5KG/9.9LBS	7KG/15.4LBS
Lifespan	100,000hrs	100,000hrs	100,000hrs
EMF Emission	-Negligible	- Negligible	- Negligible
Warranty	3 Years	3 Years	3 Years
Best For	Targeted Treatment	Targeted Treatment	Full Body Treatment
Max Area	36" x 34"	46" x 34"	52" x 34"

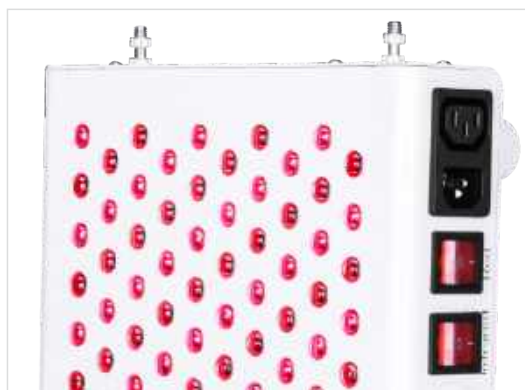


Rouge and Rouge Pro Series Irradiance

Rouge and Rouge Pro 500

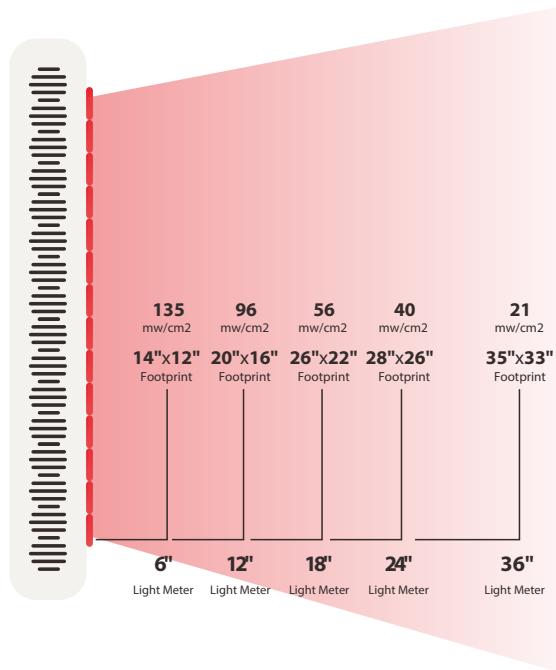


Rouge and Rouge Pro 1000

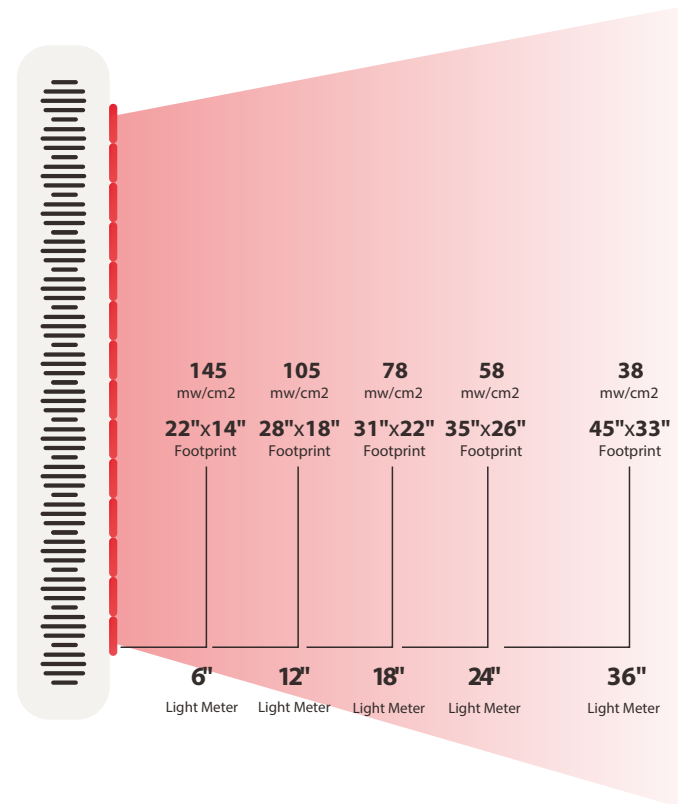

Rouge and Rouge Pro Series Irradiance

Rouge and Rouge Pro 1500

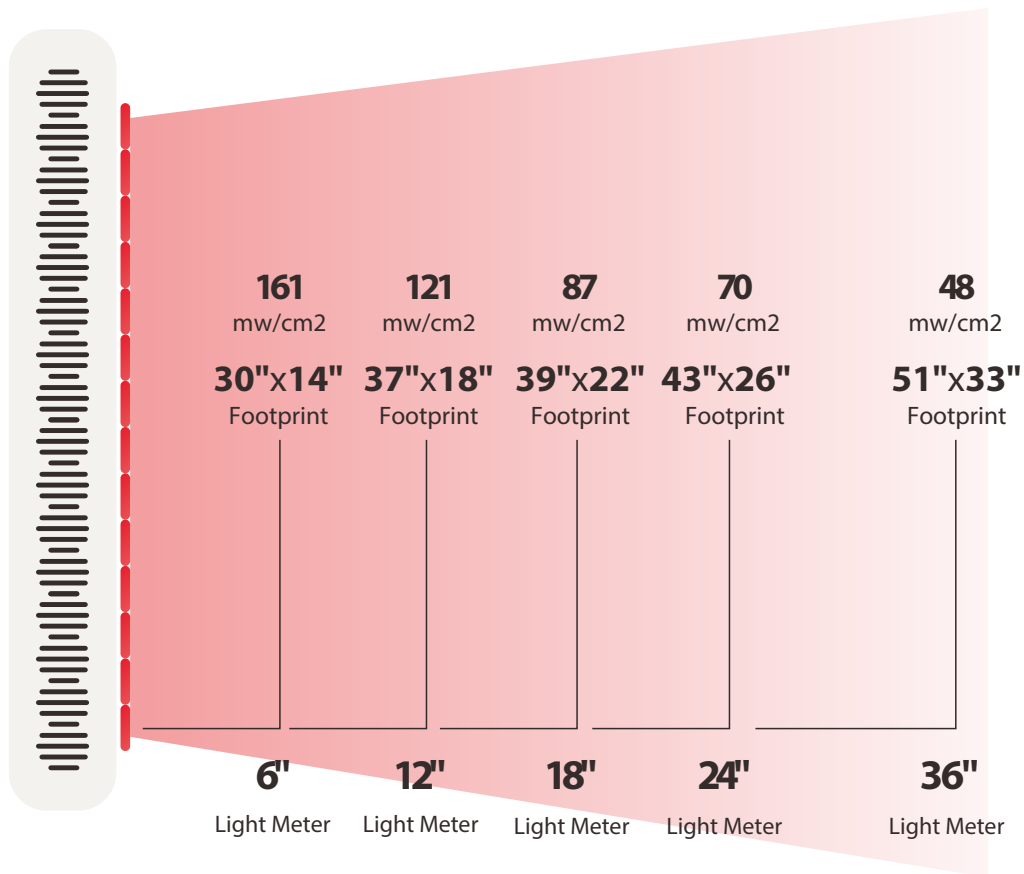


Rouge Pro Series Optional Mounting Configurations



Tabletop Stand



Door Mount

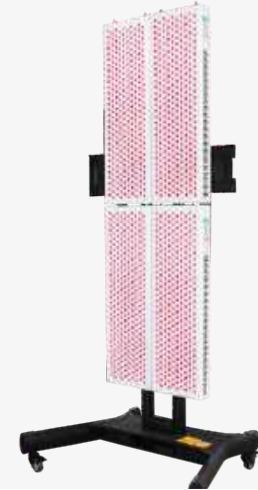

Mobile Stand

Instructions

- Distance: Apply from 6-12 inches.
- Session Time: Continue for 2-15 minutes, depending on distance.
- Frequency: Do this 2 times to 4 times a week until results are felt.
- Maintenance: Continue using several times a week.

Rouge Series Optional Mounting Configurations




Tabletop Stand



Door Mount

Instructions

- Distance: Apply from 6-12 inches.
- Session Time: Continue for 2-15 minutes, depending on distance.
- Frequency: Do this 2 times to 4 times a week until results are felt.
- Maintenance: Continue using several times a week.

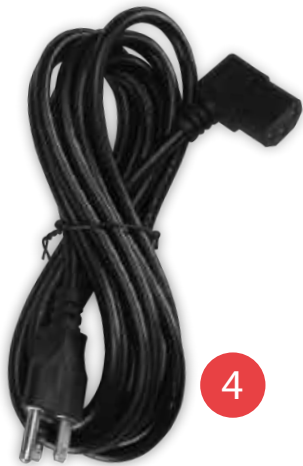

Package Contents

1 LED Therapy Light

3 Cable Assembly & Snap LinksEye

5 Protection

7 Door Hook

2 Patch Cable


4 Power Cable

6 Hanging Hardware

Install Your Therapy Light

STEP 1 ▶


Lay all of the included hardware, including your Rouge LED, on a padded surface. Next, screw in the threaded posts to the top of the device and hand-tighten with the included Allen wrench. Then, screw the cable supports to the threaded posts, as shown in PIC.1.

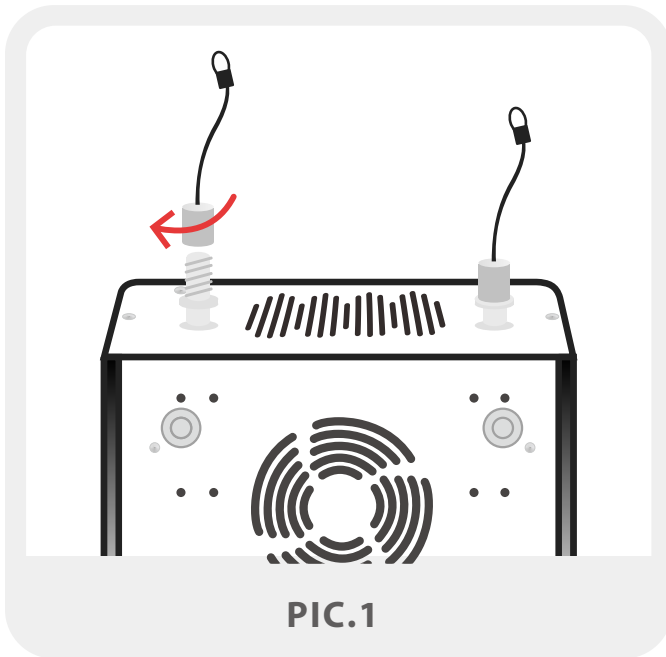


PIC.1

STEP 2 ▶

Hang the included custom hook on any standard-sized door, as shown in PIC.2.

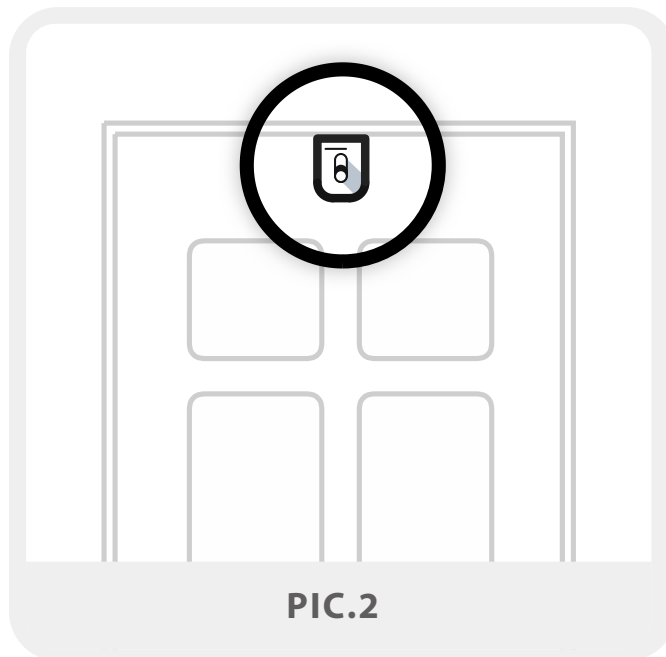


PIC.2

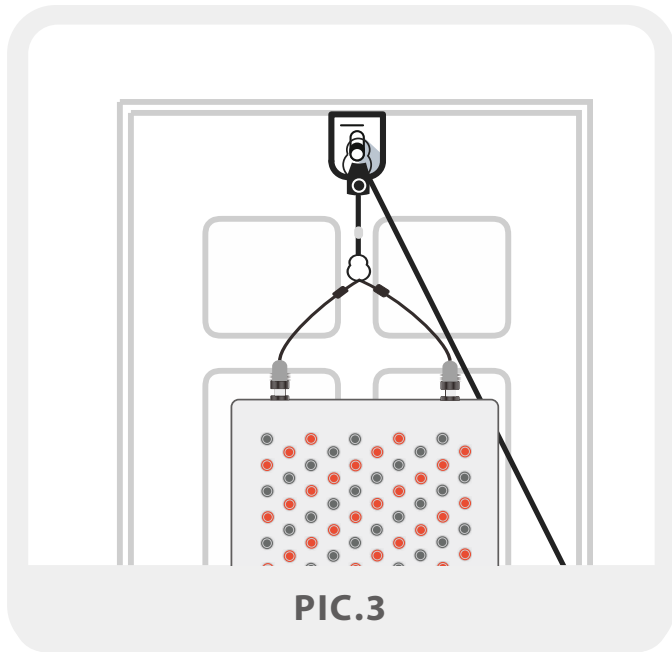

STEP 3 ▶

Connect the cable supports to the snap link attached to the height-adjusting system. Then, connect the snap link on the opposite end directly to the door hook, as shown in PIC.3.

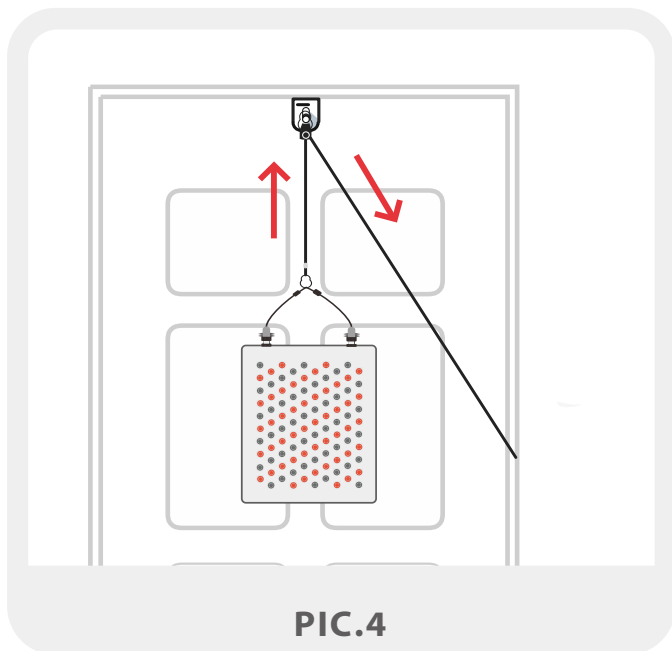


STEP 4 ▶

To raise your RougeLED to your desired height, pull down on the loose cable while supporting the device with your opposite hand, as shown in PIC.4.




Be sure to adequately support the device while raising and lowering it with the height-adjusting system .

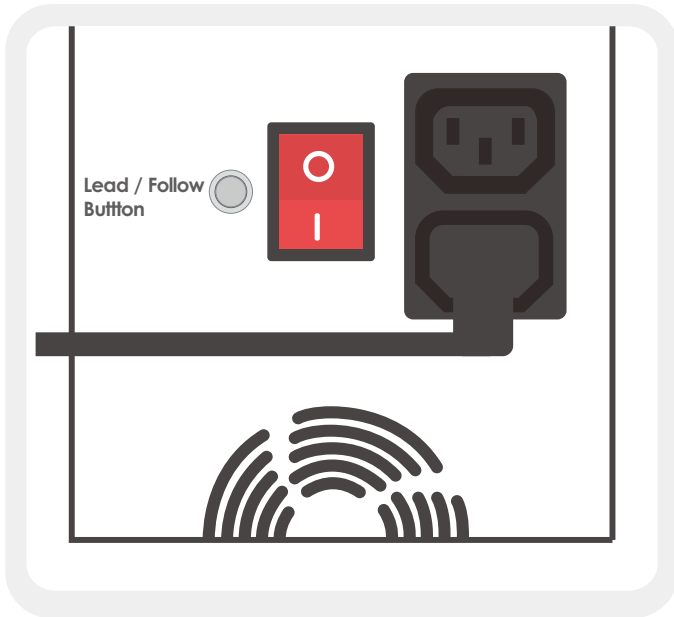

Setting up Your Rouge Pro Series Light

STEP 1 ▶

Plug the power cord in the female socket and press the red switch, it will become illuminated.

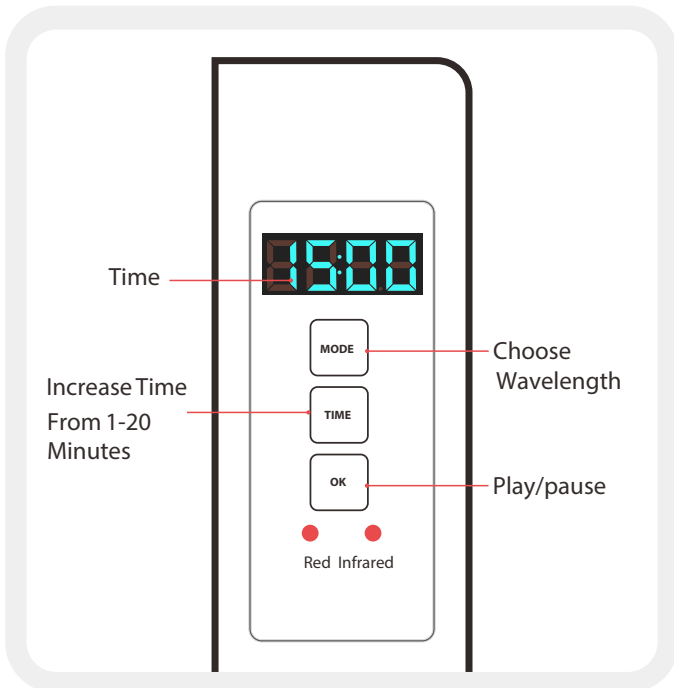


STEP 2 ▶

1. You can choose only red or only IR or both by selecting the "Mode" button. The red indicator light will illuminate indicating which wavelength configuration is selected.
2. Adjust the time by pressing the "Time" button.
3. Press "ok" to run the light, start the treatment.




Remember that near infrared LEDs emit light that is NOT visible to the naked eye.

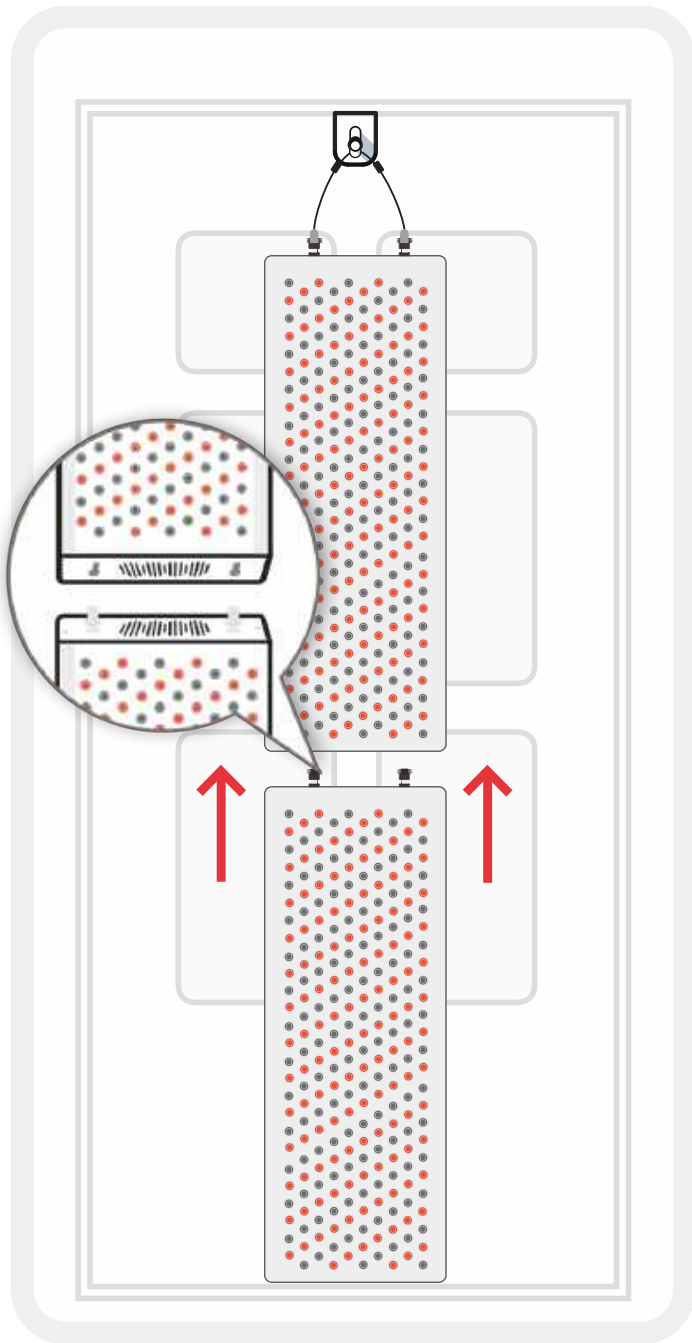

Connecting Multiple Rouge Pro Series Light Units


STEP 1 ▶

To connect the lower Rouge LED, lift it into position below the lead device that is already in place. Line up the threaded posts with the holes in the end caps of the upper device. Then, slide the lower threaded posts into the holes, and slide forward to lock in place.

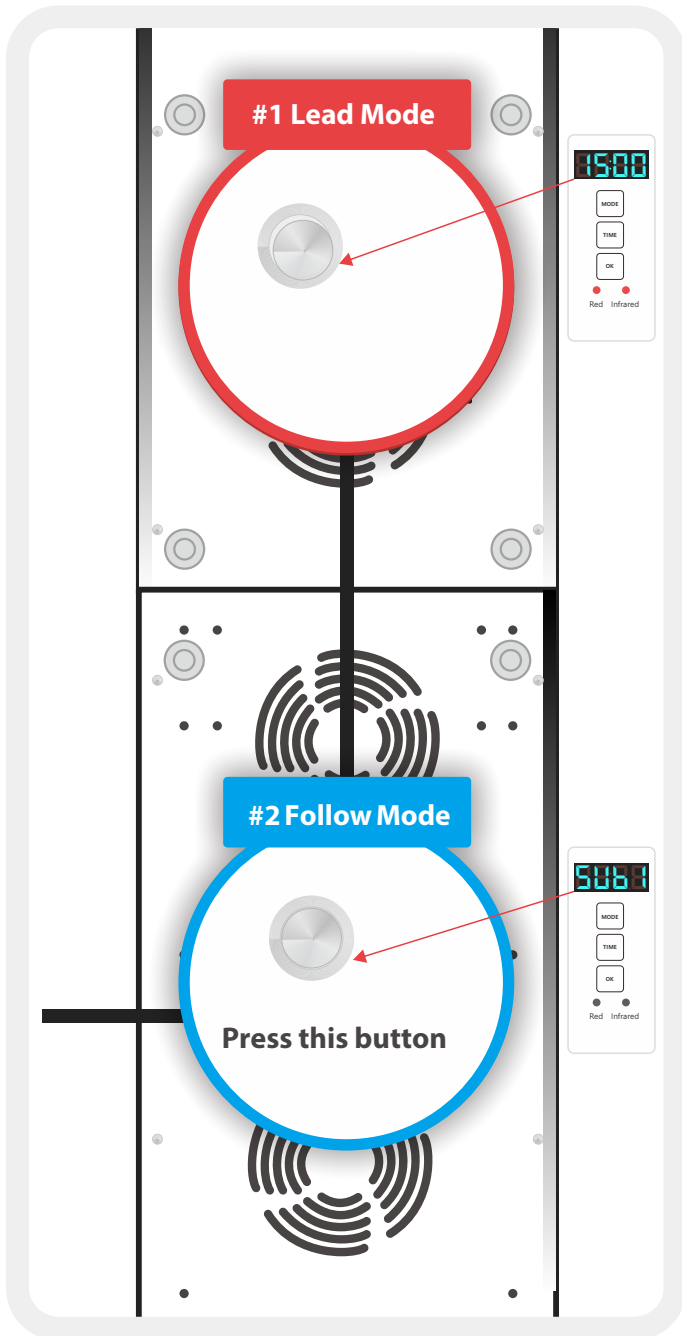

STEP 2 ▶

1. Choose #1 light as leader.
2. Press #2 light's follow button, the screen will show "SUB1", #2 lights stay in "Follow Mode".
3. You only need to control the lead light, other lights will follow the lead light's instruction.

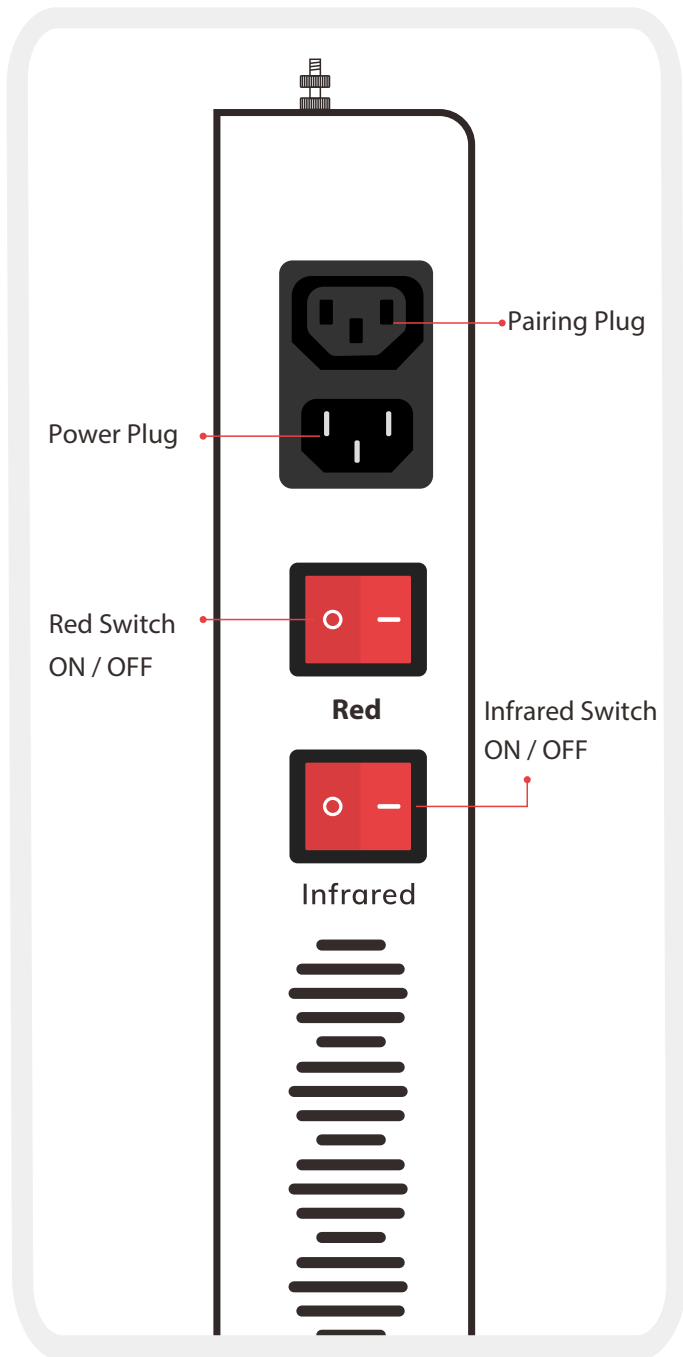


Setting up Your Rouge Series Light

STEP ▶

1. Plug the power cord into the female socket of your therapy light.
2. Plug the other end of the power cord into a wall socket.
3. Rouge Series light has two red power switches, both on the same side of the power cord socket. Each switch controls 50% of the LEDs in the light. One control all RED LEDs, and the other control Infrared LEDs.
4. Press the switches, the light fixture will work in full power.



Remember that near infrared LEDs emit light that is NOT visible to the naked eye.

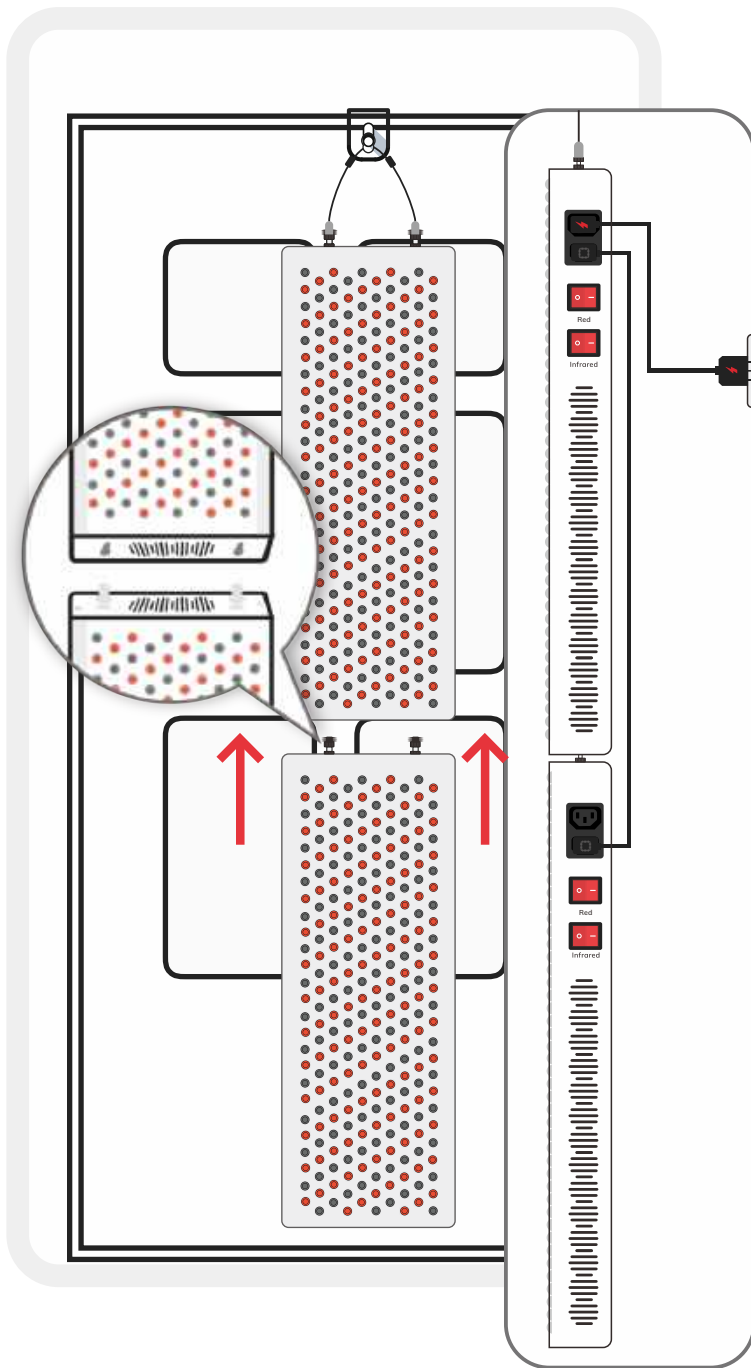

Connecting Multiple Rouge Series Light Units

STEP ▶

1. Line up the threaded posts with the holes in the end caps of the upper device. Then, slide the lower threaded posts into the holes, and slide forward to lock in place.
2. Plug the connect cable to the socket on the lights.
3. Plug the power cord into the female socket of the light, and the other end into a wall socket, then press the switches, the light fixtures will work.

Warranty Card

Contact Information

Product

Comments

Warranty Clauses

1. Rouge Red Light Therapy devices are warranted against defects in materials and workmanship for a period of THREE (3) YEARS from the purchase date.
2. Free replacements are warranted for a period of SIX (6) MONTHS.
3. Warranty does not cover the following damage(s): Modified product, improper installation or use, improper cleaning, electrical surges, or damage caused by direct exposure to natural elements like water or fire for example.
4. Warranty does not cover the loss of third party with direct or indirect damage(s)

Warning & Clean

1. Unplug fixture before cleaning.
2. Do not use cleaning sprays. Wipe with a dry or lightly damp cloth to clean.



rouge.care

✉ E-Mail: info@rouge.care

☎ Tel: 1 (833) 848-2977

📍 Rouge Red Light Therapy 2967 Dundas St. W. #16 Toronto, ON M6P 1Z2 Canada
Toronto, ON M6P 1Z2 Canada

©2020 Rouge Red Light Therapy. All rights Reserved